



Veltop Cosy 3

Instructions of mounting

Thank you for choosing the Veltop Cosy 3 model.

This product was designed, manufactured and assembled by the company CYCLO-SYSTEM in the French department of Haute-Savoie (74).


To ensure the correct use, function and long life of your new accessory, please follow the instructions for use and assembly carefully. We would like to remind you that the purpose of the cover is to protect the user from rain, cold and sun.

This version can be fitted to all wheelchairs and power wheelchairs equipped with 2 vertical tubes or a horizontal bar behind the seat.

Warning

Failure to comply with these instructions for use may impair the proper functioning of the guard and lead to the failure of parts. This would jeopardise the effectiveness of the system and could affect your safety when travelling.


Choice of 2 types of rear attachment (included)

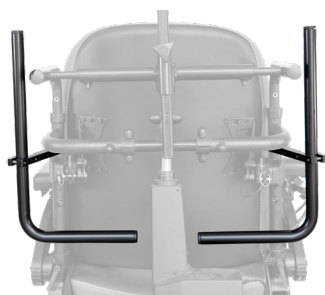


Attachment to the 2 vertical tubes



Attachment to the horizontal bar

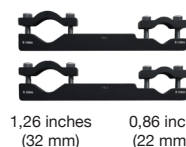

2 Brackets

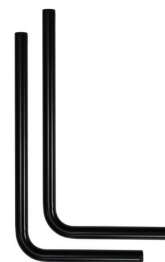

2 L-shaped tubes

Cosy 3 - A



Canopy + front part



Double open-end wrench



2 Plates



2 Alan key

Cosy 3 - B



Canopy



Double open-end wrench



2 Plates



2 Alan key



Transparent front part

Cosy 3 - C



Canopy



Double open-end wrench



2 Plates




2 Alan key



Upper part of the leg protection



Transparent front part

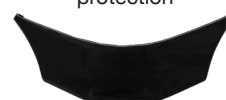

Lower part of the leg protection

Option 1 - Mounting on the 2 vertical tubes


Place the brackets (side 1,26 inches) on the two vertical tubes behind the seat and tighten them firmly.



or

Option 2 - Mounting on the horizontal bar

Place and fasten the brackets (side 1,26 inches) on the horizontal bar.

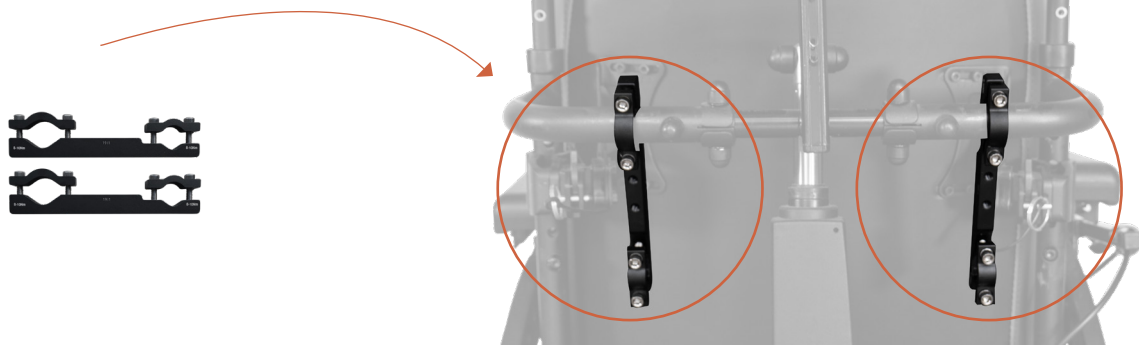

1 Fitting the L-shaped tubes

Insert the L-shaped tubes into the brackets (side 0,86 inch) and tighten the screws.



2 Fitting the plates

Insert the plates fitted with the two half-shells into the top of the tubes, ensuring that the armrests, when raised, do not interfere with the plates. The plates must be parallel to the floor and set at the same height on both sides.



3

Fitting the canopy

Slide the two plates of the canopy along the large plate and then tighten the two clamping screws.

If the seat back is permanently reclined, you can adjust the angle of the plates so that they and the hood always remain parallel to the ground.



4

Fitting the transparent front part


The transparent front part is attached using the zip fastener.



5

Fitting the upper part of the leg protection

Centre the upper part of the leg protection and attach it to the transparent part using the Velcro fastener.

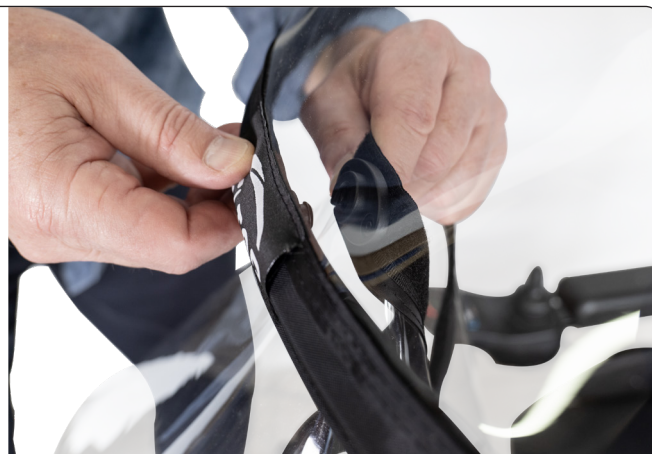

6

Fitting the lower part of the leg protection

Connect the lower part of the protection using the zip fastener to cover the underside of the legs.



Configurations



Adjustments

You can easily adjust the height of the canopy by changing the position of the large plate on the support tube.

The canopy should be as close to the head as possible.

The width between the two plates should not be greater than the width of the wheelchair.

You can also move the large plate back if it gets in the way of your arms by positioning the two half-shells at the front.



CYCLO-SYSTEM

3621 route d'Albertville
74320 SEVRIER
FRANCE

T : 0033 (0)4 50 77 56 17

F : 0033 (0)4 50 63 45 36

info@veltop.com

www.veltop.eu/en



Ref. VELCOSY3